

BAUSCH+LOMB
ULTRA[®] contact
lenses

THE FAMILY BUILT ON INNOVATION

Provide your patients all-day comfort and consistently clear vision.



The widest range of contact lenses with immediate in-office fitting*

*Based on a fit set comparison of leading brands: Biofinity, Air Optix, Acuvue Oasys, and Acuvue Vita.



Built With MoistureSeal® Technology For All-Day Comfort

💧 Maintains 95% of
lens moisture for
a full 16 hours¹



SEE
TECHNOLOGY
IN ACTION

*Results from a multisite 12-month study of Bausch + Lomb ULTRA® contact lenses when worn on a seven-day extended-wear basis. A total of 816 patients were enrolled (406 Bausch + Lomb ULTRA® and 410 control lenses) by 34 independent investigators. Performance experience ratings were obtained at follow-up visits using a 0 to 100 rating scale, with 0 denoting the most unfavorable rating and 100 denoting the most favorable rating. The results represent the average rating over all follow-up visits.


†Bausch + Lomb ULTRA® contact lenses are indicated for daily wear or extended wear of up to seven days. As with any contact lens, problems can result in serious injury to the eye, including loss of vision, and clinical studies have shown that the risk of serious adverse reactions is increased when lenses are worn overnight. It is essential that patients follow the eye care practitioner's direction and all labeling instructions for proper use of lenses and lens care products. Consult the package insert for complete information.



Advanced

MoistureSeal[®]
TECHNOLOGY

Proprietary technology uses a two-phase polymerization process¹

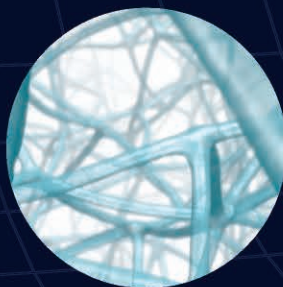

PHASE I


Silicone backbone is formed, which provides excellent oxygen transmissibility



PHASE II

PVP grows around and throughout the silicone backbone, attracting moisture

Average patient ratings during a 12-month clinical study*

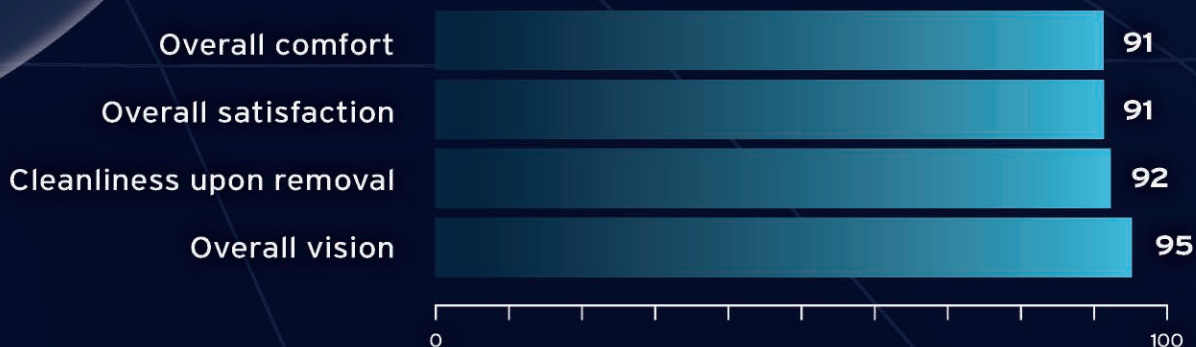

A year-long, comprehensive clinical study showed that Bausch + Lomb ULTRA[®] contact lenses provided an exceptional level of comfort when worn for extended wear.*†



Built With OpticAlign® Design For Stability

- ☉ Provides ballasting in the lower hemisphere for a stable fit*




of patients agree the lens helps reduce halos and glare even in low-light conditions†



SEE
TECHNOLOGY
IN ACTION

*Results from a 7-investigator, multisite 2-week study of Bausch + Lomb ULTRA® for Astigmatism contact lenses on 157 current soft contact lens wearers.

†Results of an online survey with patients who completed an evaluation program for Bausch + Lomb ULTRA® for Astigmatism contact lenses and wore their trial lenses for ± 4 days (n=426). Survey results include patients who strongly agreed, agreed, or slightly agreed (on a 6-point agreement scale) with the surveyed statement, with a margin of error $\pm 2.4\%$.

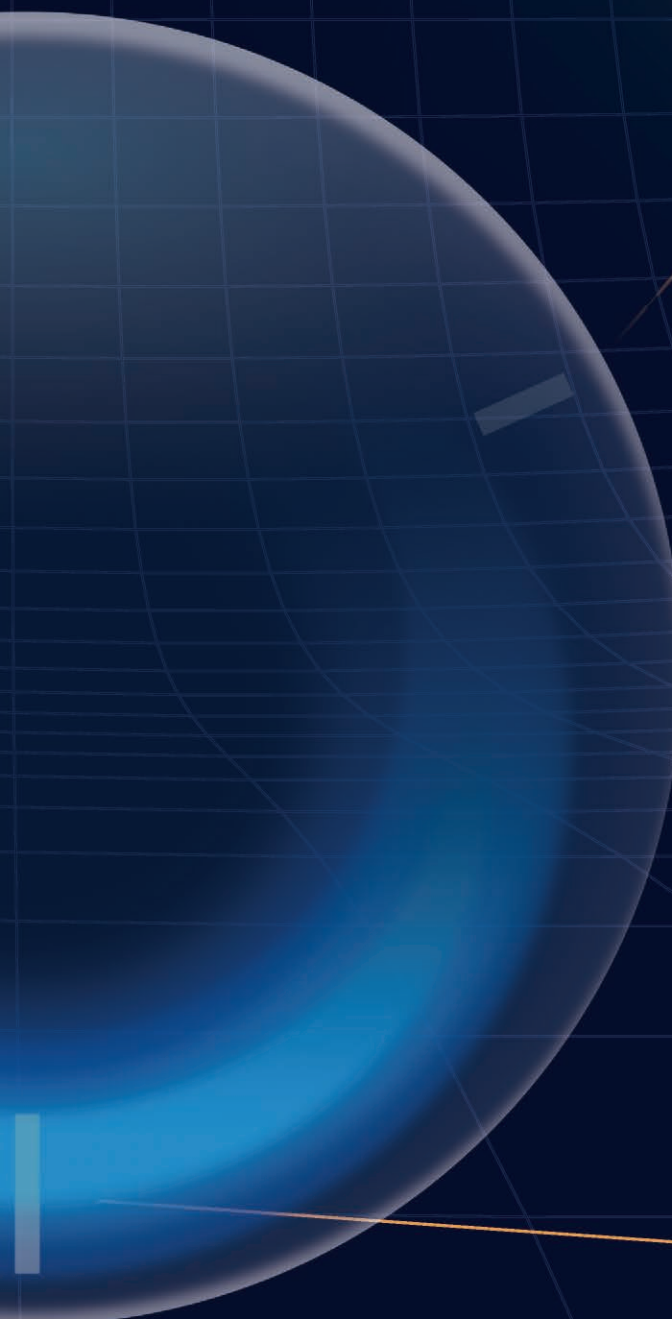

Axis indicator
designating the axis
of cylinder power

Rotation at dispensing of $\leq 5^\circ$
for 94% of patients*

**Optimized in both axes with
spherical aberration control** in 0.25D
steps to help reduce halos and glare^{†*}

**Only monthly toric with a -2.75D
Cylinder** available in a fit set
and as a standard offering

Orientation mark for
easy insertion and
measuring rotational
stability



Built With 3-Zone Progressive™ Design For Seamless Vision

- ③ Includes three distinct power zones to deliver exceptional vision and seamless transitions across all distances¹




SEE
TECHNOLOGY
IN ACTION

*When the ECP followed the fitting guide for the 3-Zone Progressive™ Design of the PureVision®2 for Presbyopia lens. Thirty-nine ECPs (from 10 countries) refitted 441 existing soft contact lens-wearing presbyopes into PureVision®2 for Presbyopia lenses. Patients returned for follow-up visits after 1-2 weeks.






of patients were successfully fit in one visit*

Optimized near and intermediate vision with the clear distance vision patients demand

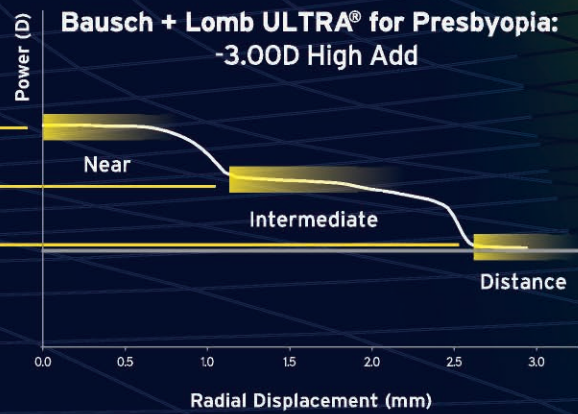

Unlike other brands, Bausch + Lomb ULTRA® for Presbyopia demonstrates consistent power in each zone of vision: near, intermediate, and distance

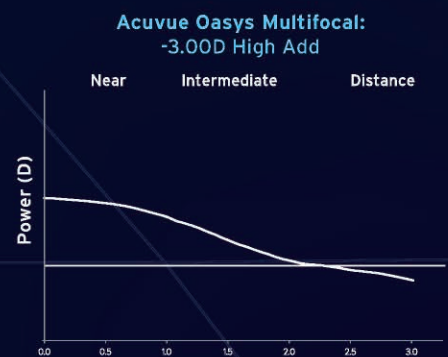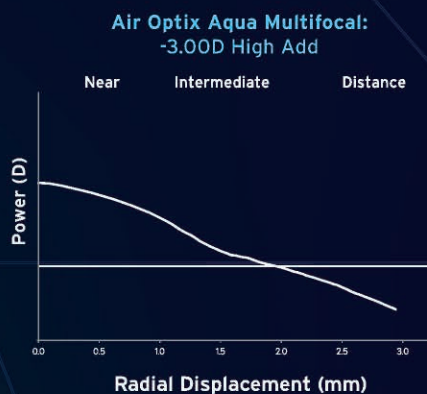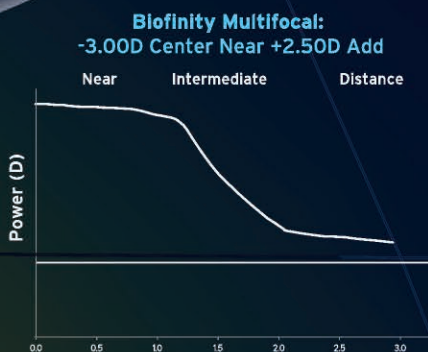



Built With Integrated Technologies For Success

All-day comfort, stable
fit, seamless vision

- 🎯 3-Zone Progressive™ Design includes three distinct power zones to deliver exceptional vision and seamless transitions across all distances¹
- 🌀 OpticAlign® Design provides ballasting in the lower hemisphere to minimize lens rotation— $\leq 5^\circ$ rotation for 95% of patients¹
- 💧 MoistureSeal® technology maintains 95% of lens moisture for a full 16 hours¹

¹Results of an online survey with Eye Care Professionals who completed an evaluation program for Bausch + Lomb ULTRA® Multifocal for Astigmatism contact lenses (n=219). Survey results include Eye Care Professionals who strongly agreed, agreed, or slightly agreed (on a 6-point agreement scale) with the surveyed statement, with a margin of error $\pm 2.8\%$.

[†]Results of an online survey with Eye Care Professionals who completed an evaluation program for Bausch + Lomb ULTRA® Multifocal for Astigmatism contact lenses (n=219). Survey results include Eye Care Professionals who reported that it took less chair time, no extra chair time, 1-2 minutes of additional chair time, or 3-5 minutes of additional chair time to fit Bausch + Lomb ULTRA® Multifocal for Astigmatism with a margin of error $\pm 5.3\%$.

Among ECPs who fit Bausch + Lomb ULTRA®
Multifocal for Astigmatism



agreed this lens allows
them to fit patients who
have previously dropped
out of contact lenses*



agreed it took no more
than 5 minutes longer to
fit this lens compared to
toric lenses†



agreed it is easy to get a
successful fit during the
first visit with this lens‡



SEE
SUCCESS
IN ACTION

‡Results of an online survey with Eye Care Professionals who completed an evaluation program for Bausch + Lomb ULTRA® Multifocal for Astigmatism contact lenses (n=219). Survey results include Eye Care Professionals who reported that successful fitting was extremely easy, easy, or somewhat easy (on a 6-point scale) during the first visit, with a margin of error ±4.1%.

Essential Savings For Better Business

ESSENTIALS

FIT YOUR LIFE

BY BAUSCH + LOMB

BUY

IN-OFFICE

ESSENTIALS



Buy 4 boxes
Bausch + Lomb ULTRA® family of lenses



HOME

ESSENTIALS HOME



Buy 4 boxes
Bausch + Lomb ULTRA® family of lenses



LITE

NEW

ESSENTIALS LITE



Buy 2 boxes
Bausch + Lomb ULTRA® family of lenses



BUNDLES

Buy 4 boxes
Bausch + Lomb ULTRA® family of lenses



Buy 2 boxes
Bausch + Lomb ULTRA® family of lenses



YOUR PATIENTS PAY ONLY FOR BAUSCH + LOMB ULTRA®.

Your Patients Benefit And Your Practice Prospers

RECEIVE FREE

Receive a 1-month trial of Biotrue® ONEday lenses + 2 x twin packs of Biotrue® multi-purpose solution.



Receive a 1-month trial of Biotrue® ONEday lenses + 2 x coupons for twin packs of Biotrue® multi-purpose solution.



Receive 1x twin pack of Biotrue® multi-purpose solution.



Receive 2x twin pack of Biotrue® multi-purpose solution.



Receive 1x twin pack of Biotrue® multi-purpose solution.



THE REST IS FREE.

*#1 Eye Care Professional Recommended Multi-purpose Solution. QuintilesIMS ProVoice Survey in Canada Nov 2017.

More Options for More Patients

				
SPHERICAL POWERS	+6.00D to -12.00D in 0.25D steps (in 0.50D steps above -6.00D)	+6.00D to -9.00D in 0.25D steps (in 0.50D steps above -6.00D)	+6.00D to -10.00D in 0.25D steps (including plano)	+4.00D to -6.00D in 0.25D steps
ADD POWERS			Low Add: +0.75D to +1.50D High Add: +1.75D to +2.50D	Low Add: +0.75D to +1.50D High Add: +1.75D to +2.50D
CYLINDER POWERS		-0.75D, -1.25D, -1.75D, -2.25D, -2.75D		-0.75D, -1.25D, -2.75D -2.25D, -2.75D
AXES:		10° to 180° in 10° increments		10° to 180° in 10° increments 10°, 20°, 70°, 80°, 90°, 100°, 110°, 160°, 170°, 180°
BASE CURVE	8.5 mm	8.6 mm	8.5 mm	8.6 mm



Recommend Biotrue® Multi-Purpose Solution—the #1 eye doctor-recommended brand*



*Based on survey data between 1/1/2021 and 12/31/2021.

REFERENCE: 1. Data on file. Bausch & Lomb Incorporated. Rochester, NY.

®/™Contact lenses are medical devices. For Professional use only.
Bausch + Lomb ULTRA® and MoistureSeal® are trademarks of Bausch & Lomb Incorporated or its affiliates.
All other product/brand names and/or logos are trademarks of their respective owners.
Bausch & Lomb Canada, Vaughan, Ontario, L4K 4B4.
© 2022 Bausch + Lomb Corporation. TP18827

BAUSCH + LOMB